

# Grouping FTICR-MS data with xcms

J. Bargsten

November 26, 2009

## Introduction

This document describes how to use *xcms* for aligning multiple MS spectra against each other.

## 1 Prerequisites

Lots of Preprocessing has to be done before the data is ready for aligning. First of all *xcms* and *MassSpecWavelet* are needed for further processing.

```
> library(xcms)
> library(MassSpecWavelet)
```

This documentation uses raw mzdata files from *msdata* as example data set. Assuming that *msdata* is installed, we locate the path of the package and extract the datafiles.

```
> library(msdata)
> mzdatapath <- system.file("fticr", package = "msdata")
> mzdatafiles <- list.files(mzdatapath, recursive = TRUE, full.names = TRUE)
> cat("Starting xcmsDirect.Rnw")
```

Starting xcmsDirect.Rnw

The *xcmsSet*-Constructor parses the given files and applies peakpicking using the *MassSpecWavelet* algorithm, leading to a *xcmsSet* object with 2 sampleclasses, ham4 and ham5, and 5 samples, respectively.

```
> data.mean <- "data.mean"
> xs <- xcmsSet(method = "MSW", files = mzdatafiles, scales = c(1,
+   4, 9), nearbyPeak = T, verbose.columns = FALSE, winSize.noise = 500,
+   SNR.method = "data.mean", snthr = 10)
```

```

HAM004_641fE_14-11-07--Exp1.extracted:
HAM004_641fE_14-11-07--Exp2.extracted:
HAM004_641fE_14-11-07--Exp3.extracted:
HAM004_641fE_14-11-07--Exp4.extracted:
HAM004_641fE_14-11-07--Exp5.extracted:
HAM005_641fE_14-11-07--Exp1.extracted:
HAM005_641fE_14-11-07--Exp2.extracted:
HAM005_641fE_14-11-07--Exp3.extracted:
HAM005_641fE_14-11-07--Exp4.extracted:
HAM005_641fE_14-11-07--Exp5.extracted:

```

## 2 Calibration

*calibrate* can be used to correct the  $m/z$  values in a `xcmsSet`. It needs a `xcmsSet` and a list of  $m/z$  value which should be found in the object. To show this on a example a sample of ham4 is created and discalibrated a bit after getting some  $m/z$ :

```

> xs4 <- xcmsSet(method = "MSW", files = mzdatafiles[1], scales = c(1,
+   4, 9), nearbyPeak = T, verbose.columns = FALSE, winSize.noise = 500,
+   SNR.method = "data.mean", snthr = 10)

```

```

HAM004_641fE_14-11-07--Exp1.extracted:

```

```

> masslist <- xs4@peaks[c(1, 4, 7), "mz"]
> xs4@peaks[, "mz"] <- xs4@peaks[, "mz"] + 1e-05 * runif(1, 0,
+   0.4) * xs4@peaks[, "mz"] + 1e-04

```

The `xcmsSet` now can be calibrated again with the  $m/z$  from the `masslist`. The plot shows the reference masses with the distances to the found ones and the regression-line.

```


> xs4c <- calibrate(xs4, wishlist = masslist, method = "edgeshift",
+   mzabs = 1e-04, mzppm = 5, neighbours = 3, plotres = TRUE)

```

```

0.0005784024
0.00060213
0.000619931

```

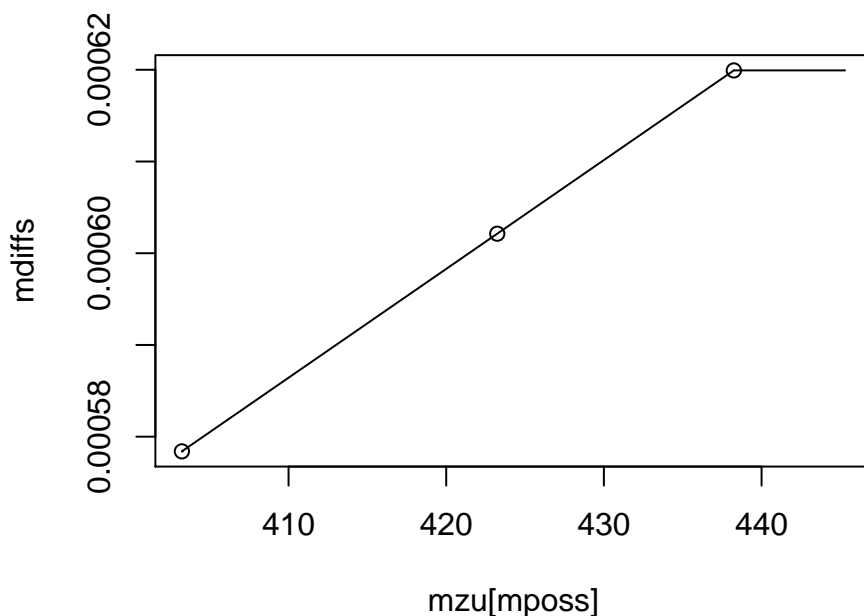

The method "shift" adds a value to each  $m/z$ , "linear" does a regression and edgshift does a regression but uses a shift before the smallest and after the biggest  $m/z$  from the wishlist.

These steps are necessary to create a usable input for *mzClust*. However, if you have already stored the data in a `xcmsSet`, you can skip the steps above.

### 3 Aligning

Now we can align `xs` with *mzClust*. The result is a clone of `xs` enhanced by the result of *mzClust*. For a description of the arguments *mzClust* takes, see helppage of the function.

```
> xsg <- group(xs, method = "mzClust")
```

```
1.14  6.82  17.05  21.59  27.27  31.82  37.50  38.64  42.05  47.73  53.41  54.55  61.36
```

```
> xsg
```

An "xcmsSet" object with 10 samples

Time range: -1--1 seconds (0-0 minutes)

Mass range: 400.1046-445.2931  $m/z$

Peaks: 88 (about 9 per sample)

Peak Groups: 18

Sample classes: ham4, ham5

Profile settings: method = bin  
step = 0.1

Memory usage: 0.0163 MB

*mzClust* stores the grouping information like the standard *group* method of *xcms* suited for retrieval via *groups* and *groupidx*. An example is shown below.

```
> groups(xsg)[1:10, ]
```

	mzmed	mzmin	mzmax	rtmed	rtmin	rtmax	npeaks	ham4	ham5
[1,]	402.2854	402.2851	402.2859	-1	-1	-1	5	0	5
[2,]	403.2365	403.2357	403.2367	-1	-1	-1	9	5	4
[3,]	405.1089	405.1087	405.1095	-1	-1	-1	4	0	4
[4,]	409.1844	409.1837	409.1845	-1	-1	-1	5	5	0
[5,]	410.1444	410.1440	410.1448	-1	-1	-1	4	0	4
[6,]	413.2672	413.2669	413.2677	-1	-1	-1	5	5	0
[7,]	423.2374	423.2363	423.2398	-1	-1	-1	3	3	0
[8,]	424.1611	424.1606	424.1615	-1	-1	-1	5	0	5
[9,]	425.1346	425.1344	425.1353	-1	-1	-1	5	0	5
[10,]	427.2681	427.2679	427.2681	-1	-1	-1	6	5	1

```
> peaks(xsg)[groupidx(xsg)[[1]]]
```

```
[1] 402.2851 402.2851 402.2851 402.2859 402.2859
```

## 4 Postprocessing

In most cases not all samples are in one group. This can be the origin of serious problems in code, which is based on e.g. *groupval*. *groupval* sets missing peaks to NA. The solution is *fillPeaks*. It changes all NA values to random noise based on the raw data file.

```
> groupval(xsg)[1, ]
```

HAM004_641fE_14-11-07--Exp1.extracted	HAM004_641fE_14-11-07--Exp2.extracted
NA	NA
HAM004_641fE_14-11-07--Exp3.extracted	HAM004_641fE_14-11-07--Exp4.extracted
NA	NA
HAM004_641fE_14-11-07--Exp5.extracted	HAM005_641fE_14-11-07--Exp1.extracted
NA	37

```

HAM005_641fE_14-11-07--Exp2.extracted HAM005_641fE_14-11-07--Exp3.extracted
                                49                                60
HAM005_641fE_14-11-07--Exp4.extracted HAM005_641fE_14-11-07--Exp5.extracted
                                70                                80

```

```
> xsgf <- fillPeaks(xsg, method = "MSW")
```

```

HAM004_641fE_14-11-07--Exp1.extracted HAM004_641fE_14-11-07--Exp2.extracted HAM004_641fE_14-11-07--Exp3.extracted

```

```
> groupval(xsgf, "medret", "into")[1:10, ]
```

```

      HAM004_641fE_14-11-07--Exp1.extracted
402.3/-1      768754.0
403.2/-1      4735257.5
405.1/-1      761632.1
409.2/-1      4158404.5
410.1/-1      726003.9
413.3/-1      6099006.3
423.2/-1      2708391.1
424.2/-1      772516.1
425.1/-1      885238.7
427.3/-1      6302089.0
      HAM004_641fE_14-11-07--Exp2.extracted
402.3/-1      1230140.4
403.2/-1      6202417.6
405.1/-1      491944.3
409.2/-1      5004546.3
410.1/-1      532868.8
413.3/-1      4950641.7
423.2/-1      1801494.2
424.2/-1      521511.4
425.1/-1      948516.6
427.3/-1      5884065.2
      HAM004_641fE_14-11-07--Exp3.extracted
402.3/-1      810120.4
403.2/-1      6117414.1
405.1/-1      650391.6
409.2/-1      4403588.2
410.1/-1      1182671.8
413.3/-1      5517709.5
423.2/-1      2826896.2
424.2/-1      347349.0
425.1/-1      376380.5

```

427.3/-1	5354053.7
HAM004_641fE_14-11-07--Exp4.extracted	
402.3/-1	568660.8
403.2/-1	5328574.1
405.1/-1	950315.1
409.2/-1	4336554.2
410.1/-1	805050.9
413.3/-1	5008541.7
423.2/-1	2427532.8
424.2/-1	481549.3
425.1/-1	961748.7
427.3/-1	5654936.8
HAM004_641fE_14-11-07--Exp5.extracted	
402.3/-1	572090.2
403.2/-1	6429028.9
405.1/-1	1452332.4
409.2/-1	4580892.8
410.1/-1	1128403.7
413.3/-1	4856606.0
423.2/-1	1856867.6
424.2/-1	499037.6
425.1/-1	662997.4
427.3/-1	5248273.8
HAM005_641fE_14-11-07--Exp1.extracted	
402.3/-1	4095293
403.2/-1	4811391
405.1/-1	2982453
409.2/-1	1196232
410.1/-1	2872023
413.3/-1	1786533
423.2/-1	1064349
424.2/-1	2995850
425.1/-1	4431535
427.3/-1	3761371
HAM005_641fE_14-11-07--Exp2.extracted	
402.3/-1	4804762.5
403.2/-1	2581183.1
405.1/-1	2268984.5
409.2/-1	1210941.2
410.1/-1	2133219.4
413.3/-1	1061103.2
423.2/-1	688353.7

424.2/-1	2556865.3
425.1/-1	3821099.0
427.3/-1	1456574.4
HAM005_641fE_14-11-07--Exp3.extracted	
402.3/-1	4657726.8
403.2/-1	2727237.5
405.1/-1	2971705.2
409.2/-1	544048.4
410.1/-1	2466625.6
413.3/-1	892797.1
423.2/-1	897205.4
424.2/-1	2567877.2
425.1/-1	4246330.7
427.3/-1	1196915.1
HAM005_641fE_14-11-07--Exp4.extracted	
402.3/-1	3755889.7
403.2/-1	2496858.9
405.1/-1	2291624.1
409.2/-1	1346778.3
410.1/-1	2980996.6
413.3/-1	982491.3
423.2/-1	1209006.2
424.2/-1	2857624.1
425.1/-1	2977003.1
427.3/-1	1227957.6
HAM005_641fE_14-11-07--Exp5.extracted	
402.3/-1	5265972
403.2/-1	2165162
405.1/-1	3009065
409.2/-1	1187547
410.1/-1	2296774
413.3/-1	1027673
423.2/-1	1136440
424.2/-1	2892810
425.1/-1	3301529
427.3/-1	1682024

The results are suited for instance for heatmaps, etc.

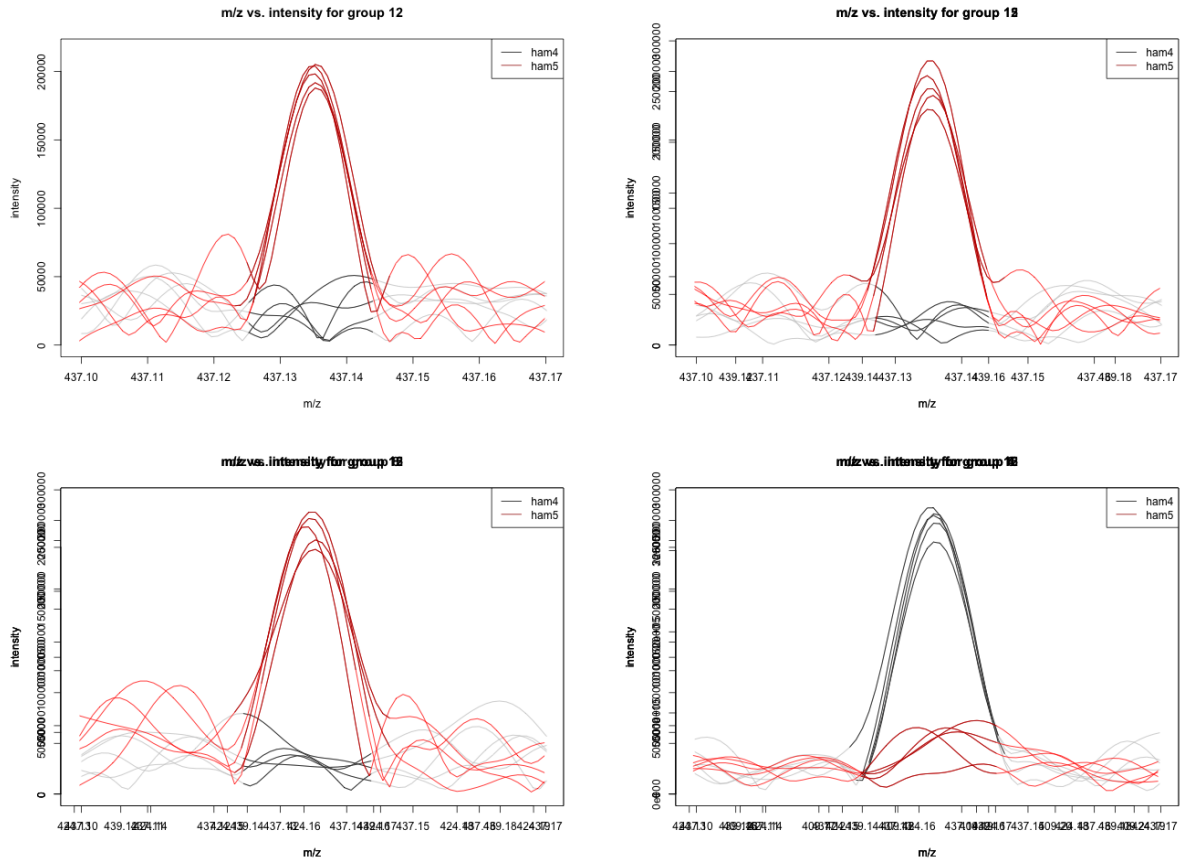


Figure 1: Auto-generated extracted spectra for the top three differentially regulated ions. Darkened lines indicate where the peaks were integrated for quantitation.

## 5 Analyzing and Visualizing Results

A report showing the most statistically significant differences in analyte intensities can be generated with the *diffreport* method. It will automatically show the superimposed peaks in the spectra for a given number of them, in this case 10. Several of those chromatograms are shown in Figure 1.

```
> reporttab <- diffreport(xsgf, "ham4", "ham5", "example", eicmax = 4,
+   h = 480, w = 640)
```

```
Processing data from sample:  1  2  3  4  5  6  7  8  9 10
```

```
group:  12 15  8  4
```

```
> reporttab[1:4, ]
```

	name	fold	tstat	pvalue	mzmed	mzmin	mzmax	rtmed
1	M437T-1_1	5.144529	20.00921	4.258336e-08	437.1353	437.1353	437.1353	-1
2	M439T-1	6.877120	25.31711	4.886837e-08	439.1508	439.1503	439.1512	-1
3	M424T-1	5.290320	19.91077	8.868895e-08	424.1611	424.1606	424.1615	-1
4	M409T-1	4.098768	-16.86272	1.555750e-07	409.1844	409.1837	409.1845	-1
	rtmin	rtmax	npeaks	ham4	ham5	HAM004_641fE_14-11-07--Exp1.extracted		
1	-1	-1	4	0	4	665915.1		
2	-1	-1	5	0	5	595855.3		
3	-1	-1	5	0	5	772516.1		
4	-1	-1	5	5	0	4158404.5		
	HAM004_641fE_14-11-07--Exp2.extracted				HAM004_641fE_14-11-07--Exp3.extracted			
1					368064.3			
2					645369.3			
3					521511.4			
4					5004546.3			
	HAM004_641fE_14-11-07--Exp4.extracted				HAM004_641fE_14-11-07--Exp5.extracted			
1					683672.7			
2					320544.9			
3					481549.3			
4					4336554.2			
	HAM005_641fE_14-11-07--Exp1.extracted				HAM005_641fE_14-11-07--Exp2.extracted			
1					2619631			
2					3586827			
3					2995850			
4					1196232			
	HAM005_641fE_14-11-07--Exp3.extracted				HAM005_641fE_14-11-07--Exp4.extracted			
1					2470892.7			
2					3606573.9			
3					2567877.2			
4					544048.4			
	HAM005_641fE_14-11-07--Exp5.extracted							
1					2683008			
2					3465884			
3					2892810			
4					1187547			